



Yes, you can.[®]

> Rea[®] Assist

To assist in maximising comfort

The Rea[®] **Assist** is ideally suited for the passive user with its comfortable seating that provides postural support. The **Assist** has been developed to facilitate the assistance that the user requires.

The improved Invacare[®] **Flex2** seat adds individual adaptability for stable and comfortable seat posture and user activity.





➤ **Innumerable possibilities in only 2 frames**

The *Rea*[®] **Assist** provides the user with an adjustable seat to ensure postural support and comfort, suitable for users seated over a long period of time. The highly adjustable seating system adjusts both in width and seat depth, with a specially designed cushion and armrest pad to provide all-round support.

➤ **Surrounding comfort**

Extra comfort is obtained by maximising contact with the chair's cushion. The deeply contoured back cushion provides excellent support to the back and sides. This specially designed Lateral cover creates optimal comfort and positioning.

➤ **Tilt-in-space for relaxation**

To optimised seating needs, **Assist** has a 20° seat tilt and an independent 30° backrest recline. These features allow the user to be active, while at the same time offering the option of easily tilting the seat if the user wishes to rest or have a short nap. The tilt and recline mechanism is operated by a gas piston or battery powered pendant control. A biomechanical backrest hinge provides support and positioning when adjusting the backrest angle.

➤ **Compatibility between product lines**

The backrest attachment has been redesigned to enhance compatibility between the *Rea*[®] **Bellis**, *Rea*[®] **Assist** and *Rea*[®] **Comfort**. The *Invacare*[®] **Flex2** range includes backrest types, cushions and other backrest options to make it easy to use the broad and numerous Invacare options for obtaining comfort.





Flex2 assortment

The interchangeable **Flex2** seat width, with optional seat solutions.



Length compensation

The leg- and backrest includes length compensation for preserving support for the head and legs through the biomechanical hinges.




Rear wheels

The large transport rear wheels roll easily, even on bumpy surfaces.



Balance

The seat has optional positions for re-balance.

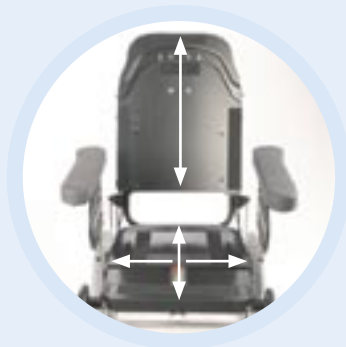

Seat size


The size of the seat is adjustable in width, depth and backrest height.



Frame length

Frame length adjustment for better stability and manoeuvring.

Technical data

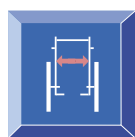


42 - 50 cm
47 - 55 cm




42 - 52 cm



47 cm

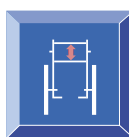

53,5 + 15 cm
(with 10 cm seat cushion)

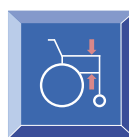


22 - 31 cm




32 - 44 cm
41 - 53 cm



0° - + 20°

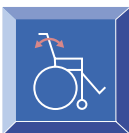

0° - 30°

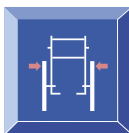


61 - 70 cm
66 - 75 cm



109 cm



101.5 cm

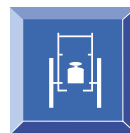

34 kg



max. 125 kg



Transport:
Weight 25 kg
73 x 65 x 48 cm



Crash tested
ISO 7176-19*

Note:

* Our wheelchairs complying with above ISO norm and with ISO 10542 and have additional features providing increased levels of occupant security and safety whilst travelling in a motor vehicle. Wheelchair users should however transfer to the vehicle seat and use the vehicle installed restraint system whenever it is feasible.

Also see:

- For compatibility of Flex2 see Flex2 brochure.
- For accessories see Accessory brochure.

Frame colour



Pearl grey

Upholstery



Grey Dartex TR23



Grey plush TR 32

All information quoted is believed to be correct at time of print.
Invacare® reserves the right to alter product specifications without prior consultation.



INVACARE® Rea AB
Box 200 · Växjövägen 303
SE-343 75 DIÖ · Sweden
Tel: +46 476 535 00 · Fax + 46 476 535 99

